virtufit

M1s
VF11001


USER MANUAL

📷 @virtufit_fitness

▶ VirtuFit

🌐 www.virtufit.com

FIG. A - MASSAGE HEAD

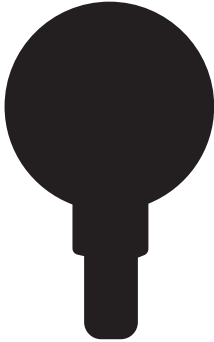Cylinder Head

EN Relax large muscle groups, such as back, legs and waist.

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Bullet Head

EN Relieves pain in soft tissue, knees and joints.

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U Head

EN Relieves stiffness of the cervical spine (neck).

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INDEX

Safety instructions	03
Guarantee	03
Use	04
Charging	04

SAFETY INSTRUCTIONS

WARNING!

Please consult your doctor before starting exercise. This is particularly important for persons with health problems. Read all instructions before using the device. VirtuFit takes no responsibility for any injuries or property damage suffered as a result of using this equipment. To reduce the risk of electric shock, fire or personal injury, this equipment should be used under the following instructions.

- For adult use only.
- Use only on a dry and clean body over clothing.
- Use only on soft tissue without causing pain or soreness. Do not use on the head or hard parts of the body or if injured.
- Press and move lightly for about 60 seconds per body part.
- Use only the attachments that give the best and desired effect.
- Bruising may occur regardless of the setting of the device or the pressure applied. Check the treated areas regularly and stop immediately in case of pain or discomfort.
- Keep fingers and hair away from the back and ventilation opening of the device, as entrapment may occur.
- Do not place objects in the ventilation opening of the device.
- Keep the device dry, free of dust and dirt and do not use water.
- Do not drop the device and use it carefully.
- Charge the device only with the original charger supplied. Do not use while charging.
- Never use or charge the appliance unattended
- Check the device and battery carefully before each use.

- Do not use the device continuously for more than one hour. Let the device rest for 30 minutes before reuse.
- Read the manual before use.

WARNING!

Do not use the device without your doctor's approval in the following cases: pregnancy, diabetes with complications such as neuropathy or retinal damage, pacemaker wear, recent surgery, epilepsy or migraines, hernias, recent joint replacements or IUDs, metal pins or plates or other concerns about your physical health. Vulnerable persons and children should be accompanied by an adult when using the device. Research is ongoing on the effects of vibration exercises for specific medical conditions. Vibration exercises in a treatment plan is advisable in some cases. This should be done on the advice of and under the supervision of a doctor, specialist or physiotherapist.

GUARANTEE

The warranty period for this item is 1 year from the day of purchase. The guarantee is only valid when this article is used for activities and normal recreational use and does not cover articles used for rental purposes.

Warranty claims are excluded if the cause of the defect is due to:

- If an item or part is deemed defective by VirtuFit, the warranty covers the repair or replacement of only the defective item or part. VirtuFit is not responsible for any costs, losses or damages incurred as a result of loss or use of this product.
- This warranty does not cover damage caused by misuse, abuse, neglect, normal wear and tear, damage caused by excessive exposure to heat or water, damage caused by improper handling and storage, damage caused by transportation, defective device caused by dropping and damage caused by anything other than defects in materials and workmanship.
- Damage exceeding the cost of the product.
- This warranty is invalid if any unauthorised repair, alteration or modification has taken place to any part of the device
- The warranty for repaired or replacement parts is valid from the date of original purchase.
- In case of improper use, neglect and/or poor maintenance.
- Failure to maintain the item in accordance with factory instructions.

USE

NOTE! The Massage Gun cannot be used while charging.

Turning on and off

- Switch the device on by pressing and holding the on/off switch on the bottom of the device for three seconds.
- The Mini Massage Gun starts in speed level 1. Press the on/off switch once to increase the speed by one level and press again to increase another level, etc.
- The maximum speed is level 4. In level 4, press the on/off switch once and the speed goes back to level 1.
- Massage the desired body part at the desired speed, applying light pressure if necessary.
- Switch off the device by pressing and holding the on/off switch for three seconds.

NOTE! To protect bones and muscles and prevent prolonged massage at the same point, the Massage Gun switches off automatically after 20 minutes of use.

Massage Head (FIG. A)

Changing the Massage Head

- Pull out and remove the message head
- Push the new message head into the device

NOTE! The massage head does not need to be rotated after pushing the message head into the device. The different message heads are explained in FIG. A.

Indicator light

- The indicator light is green when switched on
- During charging, the signal lights flash in turn. When the four signal lights are on at the same time, the battery is full.

Speed

The indicator lights represent the speed level the device is currently at.

CHARGING

NOTE! Fully charge the battery before first use. The charging input is an input for 5V-2A mobile phone USB charging cables. The device cannot be switched on or used during charging to ensure safe use. To prolong battery life, charge the device immediately when the indicator light flashes red.

Charging

- Connect the end of the supplied original USB cable to the charging port on the bottom of the device, and the other end of the USB cable can be connected to a USB adapter or USB port of a computer/laptop. The USB adapter is not included.
- When charging with a USB power supply, do not use a power supply with an output specification higher than DC 5V. Otherwise, the device will be damaged.
- Fully charge the device, which will take about 2 hours.
- After the device is fully charged for the first use, it can be used for about 5 hours. The specific usage time depends on the speed level and pressure applied.
- Disconnect the opader. The device is ready for use.

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